

Flotation device rigid

Accessoires for RollUP

RLP2070R



Emergency & Rescue Systems

Made in Austria

KOHLBRAT & BUNZ GMBH

A-5550 Radstadt

<http://www.kohlbrat-bunz.com/>

Inhalt

INTRODUCTION 2

1 SCOPE OF DELIVERY - RLP2070R 2

2 INTRODUCTION / PRODUCT DESCRIPTION 3

3 ASSEMBLY INSTRUCTIONS 4

4 WATER APPLICATION ACCESSORIES 8

Introduction

Thank you for choosing a high-quality rescue device from Kohlbrat & Bunz.

Before you use it for the first time, read these operating instructions in detail.

We expressly point out that we do not accept any liability for damage caused by improper use or disregard of our instructions / guidelines.

It may only be used in action if the user has become familiar with the correct and safe handling.

We recommend practicing the uses of this rescue device on a regular basis in order to ensure safe handling in the event of an emergency.

Important warning notices

For safety reasons, the Floating device (rigid) must be withdrawn from use immediately if:

- There are doubts about safe use, or
- The equipment has been stressed by a fall and that it may only be used again if a competent P person has given written consent;
- The device is damaged

1 Scope of delivery - RLP2070R



1 pc rigid float

2 pc fastening straps (fixed on the rigid device)

2 units (2x RLP2070R) are needed to provide the RollUP with the required buoyancy

2 Introduction / product description

These floating devices are the standard equipment for a RollUP when It is used in water.

There is one float device on the right side, and one float device on the left side of the RollUP required for the correct function.

Basically, 2 floats are the minimum requirement to create the necessary buoyancy so that the RollUP does not sink in the water.

An additional inflatable pillow increases the buoyancy a little in the water - but is not absolutely necessary.

For use in the sea (swell) a foot weight as well as a breast float is required. This "Marine Set" ensures an upright / vertical position of the RollUP in the water, and turns the RollUP automatically in the event of a possible tipping into the "head up position".

It is recommended to install the floats as soon as contact with water is possible.

**During a recovery in the water, never leave the RollUP unattended.
When the RollUP tips over - there is a risk of drowning for the patient.**

3 Assembly instructions



Fig. Symbol image for RollUP with float and pillow. (Application area inland waters)




Detailed view of how to thread the belts



place the float next to the RollUP.



In the area where it is to be mounted.

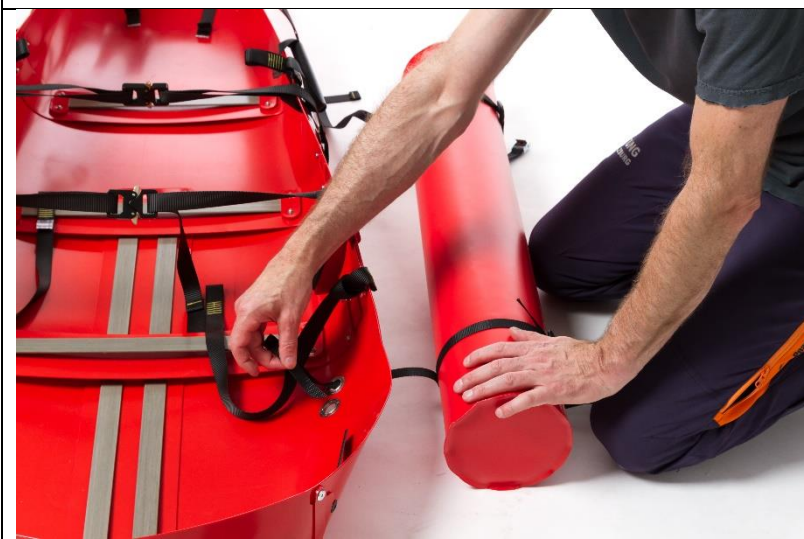

In the head area there are two eyelets. These are specially designed for mounting the float. Thread the belt through both eyelets.



pull the belt through the metal tightener and tighten it firmly

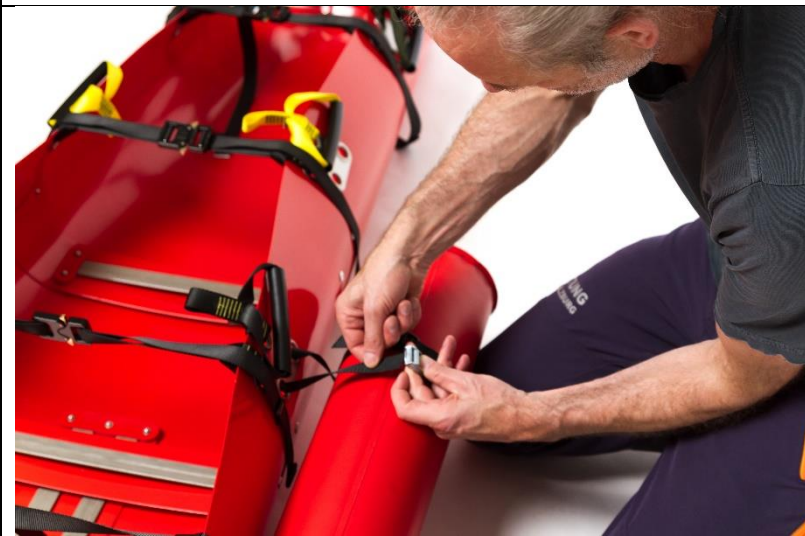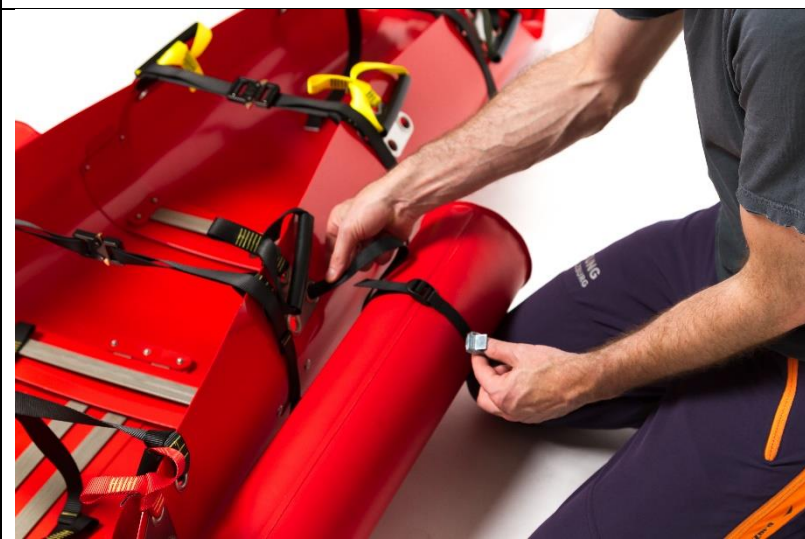

The procedure for the second belt is similar.

However, the belt is passed through two eyelets in the trapezoidal plate. It is possible that one hole is already occupied - this is not a problem simply guide the belt parallel.

Pass this strap through the metal buckle as well and tighten it.



To assemble the Pillow, (pillow Inflatable) pass the strut through the pillow and fix it.



reinsert the center brace into the RollUP



check the position and location of the pillow

4 Water application accessories

VK20X01	Pillow Pillow, Inflatable olive	
VK20X02	Pillow Pillow, Inflatable	
VK20X03	Pillow Pillow, Inflatable black cordura Easy mount	
VK20X04	Pillow Pillow, Inflatable black TPU Easy Mount	
RLP2070R	Flotation device, rigid	
RLP101	Flotation device, inflatable	
RL100G	RollUP module marine II	
RLP2077	Marine foot weight (included in the set RL100G)	
RLP102	Chest float chest inflatable (included in the set RL100G)	